



Practicon®

Gauze Guard™ Refill

70-01760 • 25/pack

- Clamps any gauze sponge over the airway for safer procedures.
- Works in tandem with all sizes of E-Prop™ Mouth Props.

Instructions (diagram on other side)

1. To ensure patient safety, tie a 24" length of dental floss into a loop and capture it within the hinge area of Gauze Guard (Figure 1).
2. Place a gauze sponge of your choice between the clamp jaws of Gauze Guard. Depending on patient's needs, gauze may be moistened and/or folded prior to loading.
3. Close clamp jaws onto gauze and press together firmly. Leading edge of gauze should align with wider end of jaws (Figure 2).
4. While holding loaded Gauze Guard in one hand and E-Prop in the other, push cylindrical portion of Gauze Guard into E-Prop as shown until it reaches the mechanical stop (Figure 3).
5. Before placing assembly into mouth, check for a secure grip of the gauze within the clamp jaws. Also check that Gauze Guard is fully inserted into E-Prop.
6. Place assembly in mouth with E-Prop over back molars and gauze positioned over airway (Figure 4). Readjust gauze as necessary to reduce gagging while maximizing airway coverage.
7. After use, dispose of Gauze Guard and gauze sponge properly. E-Prop is fully autoclavable and can be used again.

Manufactured for



PRACTICON

PRACTICAL INNOVATIONS FOR DENTISTRY

1-800-959-9505 • Fax 1-800-278-0885

practicon.com • Greenville, NC 27834

Made in China • Inner Pack 7001760_1



723896030837

Important: Gauze Guard MUST be used together with an E-Prop™ Mouth Prop (sold separately).



Do not reuse. This is a single use device.

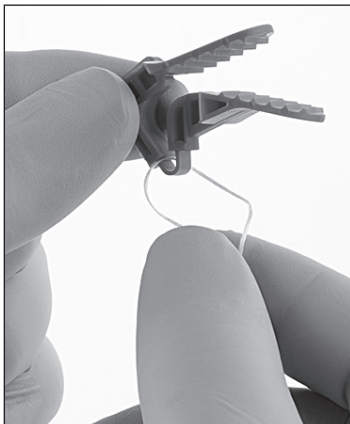


Figure 1



Figure 2



Figure 3



Figure 4